

NEW 400 HOSPIVAC

ON CASTORS

A professional surgical suction unit indicated for any surgical application guaranteeing high performances and versatility.

Thanks to its state-of-the-art maintenance-free pump it finds application in any operating theatres and general surgery.

The device is suitable for nasal, oral and endotracheal aspiration of bodily fluids (mucus or catarrh), tissues or bones from adults and children.

MAIN FEATURES



The most powerful

New Hospivac 400 is the most powerful surgical suction unit in the CA-MI line.



Non-stop operation

In order to guarantee the best performance, engine can work without interruption 24/24h.



Change-over system

Allows operator to switch to the next jar without interrupting suction.



Safety protection

3 protection systems against internal reflux.



Multi-Purpose Rail System

Flexible rail that allows the rapid change of containers.



Best Seller configuration

NEW HOSPIVAC 400 +2x 2000 mL reusable jar
in Makrolon Polycarbonate (121°C)
Product code: RE 410350

Select your configuration at page 49

The warranty lapses immediately if the unit is used without all the required accessories or with non-original CA-MI accessories.



APPLICATIONS

CE 0123

GENERAL SURGERY	NEUROSURGERY
GINECOLOGY	OBSTETRICS
OPERATING THEATRE	DENTAL PRACTICE

TECHNICAL DATA

MOTOR	Oiless and maintenance-free piston pump
POWER FEEDING	230 V / 50 Hz or 110-127 V / 60 Hz
ISO 10079-1 CLASSIFICATION	High vacuum / High flow
MAX VACUUM adjustable	-0,90 bar -90 kPa -675 mmHg value at sea level different altitudes may affect it
MAX FREE AIR FLOW RATE	90 L/min
NOISE LEVEL	46.4 dB(A)
DUTY CYCLE	Non-stop operation
WEIGHT	20 Kg
DIMENSIONS	46 x 42 x 85 cm
WARRANTY	2 years
SHIPPING CARTON	1 unit
MANUFACTURING COUNTRY	Italy



INCLUDED ACCESSORIES

Basic Setup

2 REUSABLE JARS OR DISPOSABLE SYSTEM
depending on configuration

1 SET OF AUTOCLAVABLE SILICONE TUBES
ø 8/14 mm

1 ANTIBACTERIAL & HYDROPHOBIC FILTER
single-patient

1 SAFETY TRAP BOTTLE
220 mL

1 BI-CONICAL CONNECTOR
ø 10/11/12 mm

2 JAR ACCOMODATION RINGS
size depending on the model

1 AIR SUCTION INLET

1 POWER CORD WITH SCHUKO PLUG

FS Setup

ALL THE ABOVE-MENTIONED ACCESSORIES

+ 1 FOOT-SWITCH CONTROL
with intermittent or continuous operation

Full Setup

ALL THE ABOVE-MENTIONED ACCESSORIES

+ CHANGE-OVER SYSTEM
from jar to jar by soft-touch keys

+ 1 ANTIBACTERIAL & HYDROPHOBIC FILTER
single-patient

+ 1 SET OF AUTOCLAVABLE SILICONE TUBES
ø 8/14 mm

+ 1 SAFETY TRAP BOTTLE
220 mL

+ 1 BI-CONICAL CONNECTOR
ø 10/11/12 mm

+ 1 AIR SUCTION INLET

OPTIONAL ACCESSORIES

Reusable Jars in Makrolon®, APEC® Polycarbonate or Polysulfone

AVAILABLE SIZES:
2000 mL / 4000 mL / 5000 mL

AVAILABLE SETS:

- › Complete set (Jar+cover with overflow system+ O-ring)
- › Jar (without cover)
- › Cover with overflow system
- › O-ring

Further details at page: 41



OPERATING WITH:

- REUSABLE COMPLETE TUBES' SET
- 3 silicone tubes (ø 8/14mm, L: 150cm/ 18cm each)
- + 1 bi-conical connector ø 10-11-12mm
- + 1 Antibacterial and hydrophobic filter




Disposable Flovac® System

AVAILABLE SIZES:
2000 mL / 3000 mL

AVAILABLE SETS:

- › Reusable Flovac® Canister
- › Disposable Flovac® Liner
- › Male Connectors for Flovac® Liner (5's pack)

Further details at page: 43

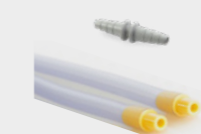

OPERATING WITH:

- FLOVAC COMPLETE TUBES' SET
- 2 tubes (ø 8 / 14mm, L: 140cm / L: 26cm)
- + 1 bi-conical connector (ø 10/11/12mm)



Rolls of silicone tube

AVAILABLE SIZES:

- › ø 8 / 14mm, Length = 1m
- › ø 8 / 14mm, Length = 10m
- › ø 8 / 14mm, Length = 50m



Cannulas

AVAILABLE MODELS

- › Yankauer Handle Flat Tip With Hole
- › Yankauer Handle Crown Tip With Hole
- › Yankauer Tube L= 180cm

Foot-switch control

with plug-in system



Multi-purpose rail & accessories

Further details at page: 12

